



FINAL REPORT

BACTERIAL FILTRATION EFFICIENCY  
AND DIFFERENTIAL PRESSURE

PROCEDURE NO. STP0004 REV 02

LABORATORY NO. 475749

PREPARED FOR:

HAN BING DOU  
SHANGHAI DOCHEM INDUSTRIES CO., LTD.  
SUITE 705, NO. 7, FU LI BUILDING, LANE 1306 JIANG NING ROAD  
SHANGHAI 200060  
CHINA

SUBMITTED BY:

NELSON LABORATORIES, INC.  
6280 S. REDWOOD RD.  
SALT LAKE CITY UT 84123-6600  
801-290-7500

Page 1 of 6

**BACTERIAL FILTRATION EFFICIENCY  
AND DIFFERENTIAL PRESSURE**

LABORATORY NUMBER:	475749
PROCEDURE NUMBER:	STP0004 REV 02
SAMPLE SOURCE:	Shanghai Dochem Industries Co., Ltd.
SAMPLE IDENTIFICATION:	Refer to Table 1
DEVIATIONS:	None
SAMPLE RECEIVED DATE:	15 May 2009
LAB PHASE START DATE:	18 May 2009
LAB PHASE COMPLETION DATE:	27 May 2009
REPORT ISSUE DATE:	28 May 2009

**INTRODUCTION:**

This test procedure was performed to determine the bacterial filtration efficiency (BFE) of various filtration materials, employing a ratio of the bacterial challenge counts to sample effluent counts, to determine percent bacterial filtration efficiency (%BFE). This procedure provides a more severe challenge to most filtration materials than would be expected in normal use. This test procedure allowed a reproducible bacterial challenge to be delivered to test materials. This method complies with ASTM F2101.

The differential pressure ( $\Delta P$  or Delta P) test determined the air exchange differential of the porous materials. The technique involved a simple application of a basic physical principle employing a water manometer differential upstream and downstream of the test material, at a constant flow rate. A digital manometer may be used in place of a water manometer.

**ACCEPTANCE CRITERIA:**

The BFE control average must be  $2200 \pm 500$  colony forming units (CFU). A BFE run with a control average of less than 1700 shall be unacceptable. Challenges greater than 2700, but less than 3000, are, in our experience, valid. Acceptance of runs with control averages exceeding 2700 shall be at the sponsor's approval.

The mean particle size (MPS) of the challenge aerosol must be maintained at  $3.0 \pm 0.3 \mu\text{m}$ .

The average % BFE for the reference material must be within the upper and lower control limits established for the BFE test.

The average Delta P result for the reference material must be within the upper and lower control limits established for the Delta P test.

#### SAMPLE PREPARATION:

BFE test samples were conditioned for a minimum of 4 hours at  $21 \pm 5^{\circ}\text{C}$  and  $85 \pm 5\%$  relative humidity prior to testing.

#### TEST PROCEDURE:

A culture of *Staphylococcus aureus*, ATCC #6538, was diluted in 1.5% peptone water (PEPW) to a precise concentration to yield challenge level counts of  $2200 \pm 500$  CFU per test sample. The bacterial culture suspension was pumped through a 'Chicago' nebulizer at a controlled flow rate and fixed air pressure. The constant challenge delivery, at a fixed air pressure, formed aerosol droplets with a MPS of approximately  $3.0 \mu\text{m}$ . The aerosol droplets were generated in a glass aerosol chamber and drawn through a six-stage, viable particle, Andersen sampler for collection. The collection flow rate through the test sample and Andersen sampler was maintained at 28.3 Liters per minute (Lpm) (1 cubic foot per minute (CFM)). Test samples, positive controls and reference material received a one minute challenge followed by a one minute vacuum cycle.

The delivery rate of the challenge also produced a consistent challenge level of  $2200 \pm 500$  CFU on the test control plates. A test control (no filter medium in the airstream) and reference material are included after 5-11 test samples. The Andersen sampler, a sieve sampler, impinged the aerosol droplets onto six soybean casein digest agar (SCDA) plates based on the size of each droplet. The agar plates were incubated at  $37 \pm 2^{\circ}\text{C}$  for  $48 \pm 4$  hours and the colonies formed by each bacteria laden aerosol droplet were counted and converted to probable hit values using the hole conversion chart provided by Andersen. These converted counts were used to determine the average challenge level delivered to the test samples. The distribution ratio of colonies for each of the six agar plates were used to calculate the MPS of the challenge aerosol.

The  $\Delta\text{P}$  test simply measured the differential air pressure on either side of the test sample using an incline, "U" tube, or digital manometer. Testing was conducted at a flow rate of 8 Lpm (volumetric).

#### RESULTS:

The results are summarized in Table 1.

Shanghai Dochem Industries Co., Ltd.  
Lab Number 475749

Bacterial Filtration Efficiency  
and Differential Pressure

The filtration efficiencies were calculated as a percent difference between test sample runs and the control average using the following equation:

$$\% BFE = \frac{C - T}{C} \times 100$$

Where: C = Average of control values.  
T = Count total for test material.

The  $\Delta P$  values were reported in mm water/cm<sup>2</sup> of test area and calculated using the following equation:

$$\Delta P (\Delta P) = \frac{\bar{M}}{\text{Test Area}}$$


Where:  $\bar{M}$  = Average mm water of test replicates.

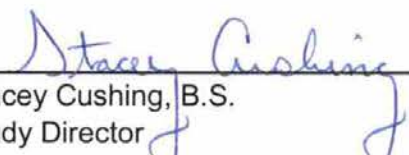
The sample holder used in the  $\Delta P$  test has a test area of 4.9 cm<sup>2</sup>.

At least one reference material is included with each set of test samples. The  $\Delta P$  values for the reference material are also recorded on control charts. The individual  $\Delta P$  values must be within the upper and lower control limits ( $\pm 3$  standard deviations) for the test. Testing met the acceptance criteria previously stated in this report.

#### STATEMENT OF UNCERTAINTY:

If applicable, the statement of uncertainty is available to sponsors upon request.

  
\_\_\_\_\_  
Technical Reviewer

  
\_\_\_\_\_  
Stacey Cushing, B.S.  
Study Director

  
\_\_\_\_\_  
28 May 2009  
Study Completion Date

anh

TABLE 1. Results  
Sample Identification: flat face mask, #1-#7

UNIT NUMBER	$\Delta P$ (mm H <sub>2</sub> O/cm <sup>2</sup> )	PERCENT BFE
1	2.6	99.8%
2	2.8	99.8%
3	2.7	99.9%
4	2.6	99.9%
5	2.5	99.8%
6	2.4	99.8%
7	2.6	99.8%

CONTROL AVERAGE: 1770 CFU

MEAN PARTICLE SIZE: 3.1  $\mu\text{m}$



Shanghai Dochem Industries Co., Ltd.  
Lab Number 475749

Bacterial Filtration Efficiency  
and Differential Pressure

All reports and letters issued by Nelson Laboratories, Inc. are for the exclusive use of the sponsor to whom they are addressed. These results relate only to the samples tested. Reports may not be reproduced except in their entirety. No quotations from reports or use of the corporate name is permitted except as expressly authorized by Nelson Laboratories, Inc. in writing. The significance of any data is subject to the adequacy and representative character of the samples tendered for testing. Nelson Laboratories, Inc. warrants that all tests are performed in accordance with established laboratory procedures and standards. Nelson Laboratories, Inc. makes no other warranties of any kind, express or implied. Nelson Laboratories, Inc. expressly states that it makes no representation or warranty regarding the adequacy of the samples tendered for testing for any specific use of application, that determination being the sole responsibility of the sponsor. Nelson Laboratories' liability for any loss or damage resulting from its actions or failure to act shall not exceed the cost of tests performed, and it shall not be liable for any incidental or consequential damages. Nelson Laboratories, Inc. has a ten year record retention policy.